

What's New from Best Medical?

Best® Total Solutions™ for Breast Cancer Diagnosis & Treatment



**Best® Breast
Localization
Needles**

- Best localization needles are available in 5, 7, 10, 15 and 20 cm lengths
- Available in all ranges of seed activities (0.08 to 0.4 mCi) with a shelf life of 90 days
- Needles are supplied sterile and are individually packaged and shielded

Best® HDR Afterloader *

- Provides 18 channels for dose delivery
- Battery back-up in case of power failure and automated wire recovery
- Source wire tracking using sensors and encoders
- Treatment data maintained in case of failure and treatment resumed after fix
- Quick source replacement process thereby reducing down time
- Available with Cobalt-60 or Iridium-192 sources. Also, Ytterbium-169 source available soon!



** Pending regulatory approval for sale in USA and other countries.*

Best® Dual-Balloon Breast Brachytherapy Applicator



- Improved dose distribution and conformity
- Less dose to critical tissues such as skin, lung, heart, chest wall, etc.
- Configured for optimal dose delivery
- Convenient to use
- US FDA 510(k) pre-market approved



© 2015 Best Medical International, Inc.

Best Medical International, Inc. 7643 Fullerton Road, Springfield, VA 22153 USA
703 451 2378 800 336 4970 www.bestmedical.com www.teambest.com

AFRICA | ASIA | EUROPE | LATIN AMERICA | MIDDLE EAST | NORTH AMERICA




Best® Total Solutions™ for Breast Cancer Diagnosis & Treatment

Best® Sonalis™ Ultrasound Imaging System

The Best® Sonalis™ Ultrasound Imaging System provides superior visualization of HDR, LDR, RF or Cryosurgical procedures.

- Totally sealed, self-healing antimicrobial keyboard with SensoFoil® Technology
- High resolution breast probe provides for needle aspiration, marker insertion and HDR applicator placement
- Longitudinal array on rectal probe provides for 140 mm length of view
- Simultaneous imaging of transverse and sagittal planes
- PC Based System provides a platform for future upgrades and application-specific modules

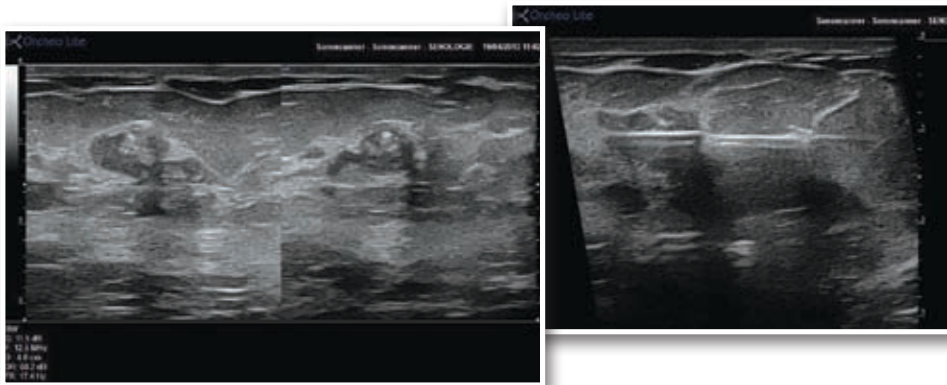

Best® Treatment Planning System (TPS)



Best® TPS is the planning system of choice for enhanced patient outcomes.

- Can be used for treatment planning with the breast balloon applicator
- Dual Activity Module creates treatment plans using multiple activity sources
- Slice Shifter™ provides ability to correct image slice positioning errors
 - Concurrent 2D & 3D Visualization allows instant visualization of seed placement, dose distribution & anatomical structures
- Pattern Load Feature – User-definable seed patterns can be used to initiate new plans

Best® UCT-1000 Universal Chair


- 4-Function Hand-Held Pendant (Control)
- Electronic variable height control/back support/foot support & horizontal tilt 12 degrees – head or foot
- Padded variable height intravenous arm rest
- Removable Stainless Steel IV Pole
- Stainless steel bilateral cushioned knee supports



Best® Brachytherapy Kit

- Catheter tubes are available in 4 colors (purple, green, clear or blue) with either radiopaque or clear nylon buttons
- Half Moon Buttons are available in radiopaque (pictured) or clear nylon

5 Single
Leader
Catheters

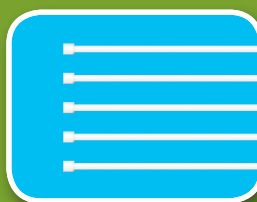

5 Half Moon Buttons



5 Friction Cuffs



5 Red Caps



5 Implant Needles

5 Stylets with Hubs



Flexi Needle Tray

Best® Flexi Needles



15G Flexi
Needle



16G Flexi
Needle



17G Flexi
Needle



Stainless Steel or Tungsten Obturator



Flexi Needle with Female Hub



BATCAM™ Multiprobe Image-Guided Radiation Therapy

BATCAM™ uses the latest compact ultrasound technology to provide increased sensitivity and resolution which yields more clinical information and expanded applications. Plus, it provides the ability to orient the probe in any rotation or translation position.

- Imports treatment plans from all major treatment planning systems
- Enhanced positioning precision using triangulated positioning techniques
- ImageSync™ Interactive Alignment provides improved visualization
- Position verification with image capture ensures accuracy
- User-friendly software components featuring a touch-screen interface




Best[®] medical *international*

Best[®] Total Solutions[™] for Breast Cancer Diagnosis & Treatment


- We manufacture both I-125 and Pd-103 seeds in the U.S. and supply directly to hospitals
- Best Medical is the only company that makes custom seeds and strands to your exact specifications — any source, any configuration, any time!
- Best[®] seeds have the best visualization in ultrasound, fluoro, X-ray, CT and MRI, facilitating real-time dosimetry
- We help set up new programs with licensing, procurement of necessary equipment, proctoring, reimbursement consulting, patient education and marketing



Best[®] Iodine-125 Seed

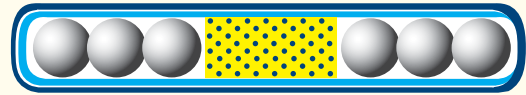

Best[®] Palladium-103 Seed

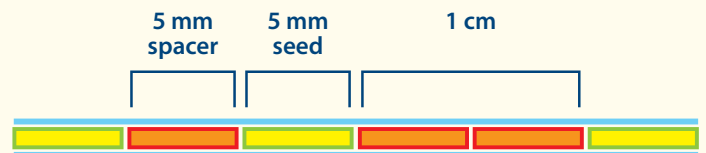

Do You Know? Best Medical is also the only supplier of Iridium-192 in nylon ribbons and Gold-198 seeds!

24 hours

We ship your order
within 24 hours, 7 days a week
Anywhere in the world!

Best[®] Fiducial Marker Kit



Stranded

Loose

Best Medical can provide any standard or customized marker configuration made to your exact specifications, sterile or non-sterile, shipped within 24 hours!

www.BestTotalSolutions.com

Best Medical International, Inc. 7643 Fullerton Road, Springfield, VA 22153 USA
703 451 2378 800 336 4970 www.bestmedical.com www.teambest.com

AFRICA | ASIA | EUROPE | LATIN AMERICA | MIDDLE EAST | NORTH AMERICA

