




Best™ Brachytherapy Kit



5 Implant Needles

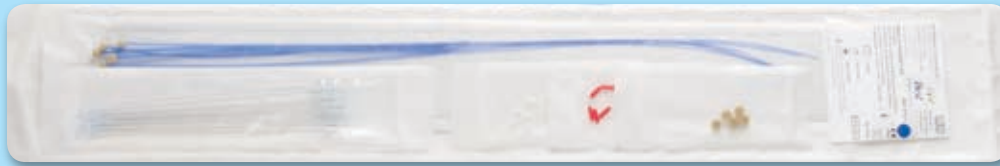

5 Stylets with Hubs




5 Friction Cuffs




5 Half Moon Buttons*



5 Red Caps

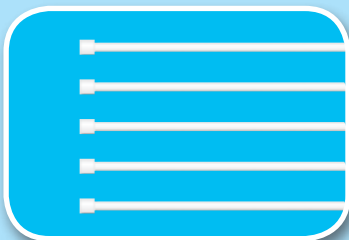

5 Single Leader Catheters**

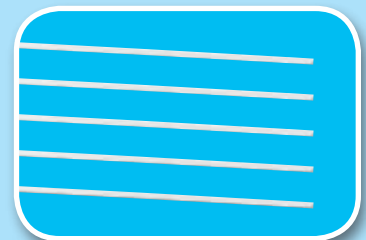

* Half Moon Buttons are available in radiopaque (pictured) or clear nylon.

** Catheter tubes are available in 5 colors (clear, blue, green or purple) with either radiopaque or clear nylon buttons.

Best™ Flexi Needles




Flexi Needle Tray

Flexi Needles are sharp or blunt, with closed-ended needle guides, and have a nominal 0.3 cm long tip. They are made of biologically compatible material. (Custom lengths are available upon request.)

Intended Use:

Flexi Needles are inserted transcutaneously, using special obturators. They are designed to be afterloaded with Iridium-192 or Iodine-125 ribbon sources. Two types of caps are available (Luer Cap or Slip Cap) to place over the loaded Flexi Needle to hold the ribbons in place during irradiation. Simulator sources can be used for localization; obturators are left in place until the time each is afterloaded with the appropriate radioactive source.

Best™ Template Needles (Flexi & Stainless)

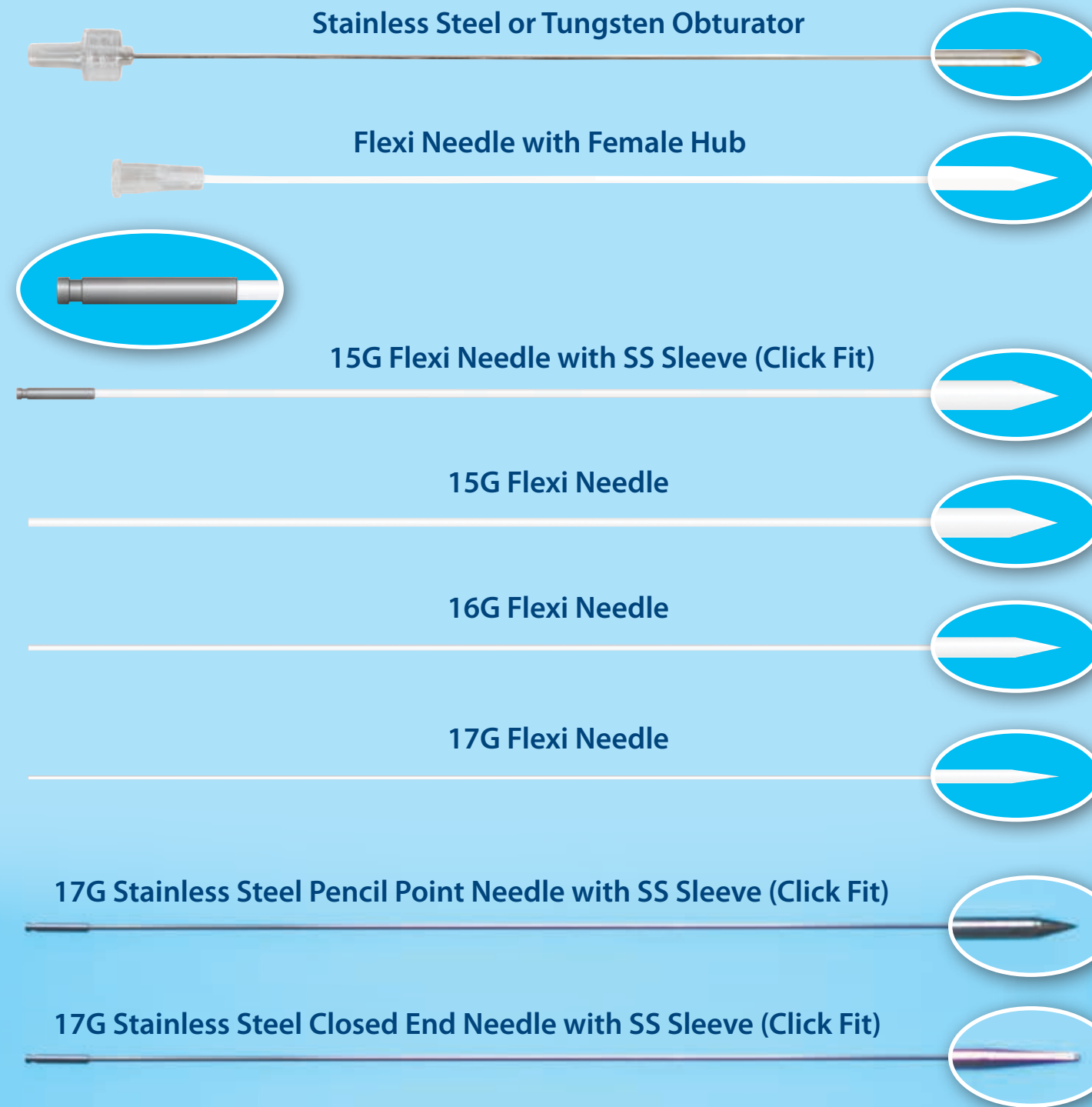


Best™ HDR Needles

These 17 Gauge HDR Needles are used with the templates, which improve the homogeneity of radiation. Needles of different lengths, closed end or open end, with or without collar, with or without funnel, are available.



Implant Needles are made of Stainless Steel and are hollow with different types of bevel at one end. Various lengths in different gauges are available, though the most commonly used needle is the 17G. Custom lengths are available upon request.

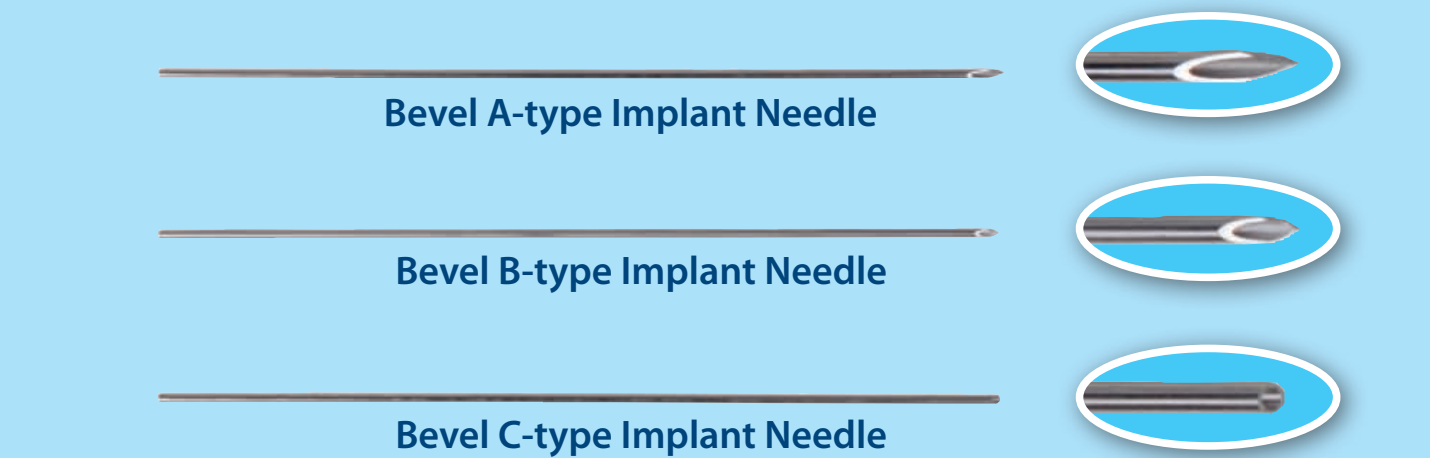


Anderson Marker Kit

The Anderson Marker Kit contains a 25 cm long implanter designed to limit implantation to 1 cm and stylet.



Adjustable Marker Kit

The Adjustable Marker Kit contains a 20 cm long implanter with adjustable stop, stylet, one spare needle and one Allen wrench.

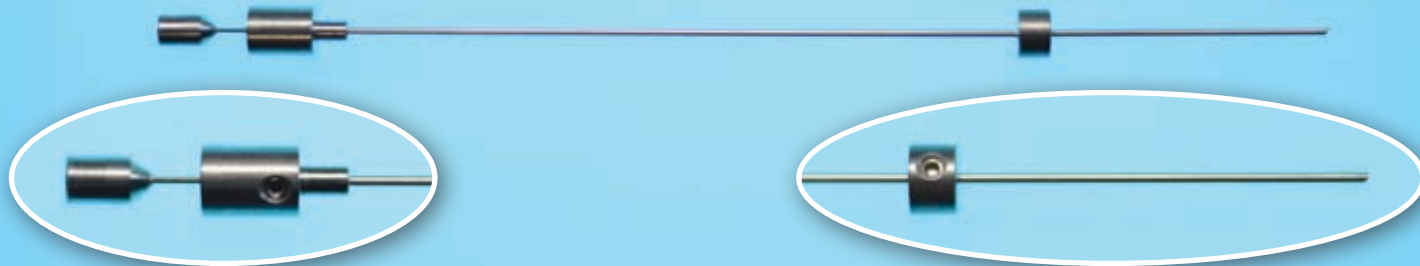


Stainless Steel Markers are sold separately. The markers are easily identified on visualization films and also aid in treatment planning. Gold Markers are also available to use with the Anderson and Adjustable Marker Kits.

Best™ Short/Long Curved Needles

Curved Stainless Steel Implant Needles have an outer diameter of 0.095 inches and inner diameter of 0.077 inches.


Two lengths are available: 3.5 inches and 7.25 inches. Each handle is 4 inches long.

Also available: Thin wall curved short needles.

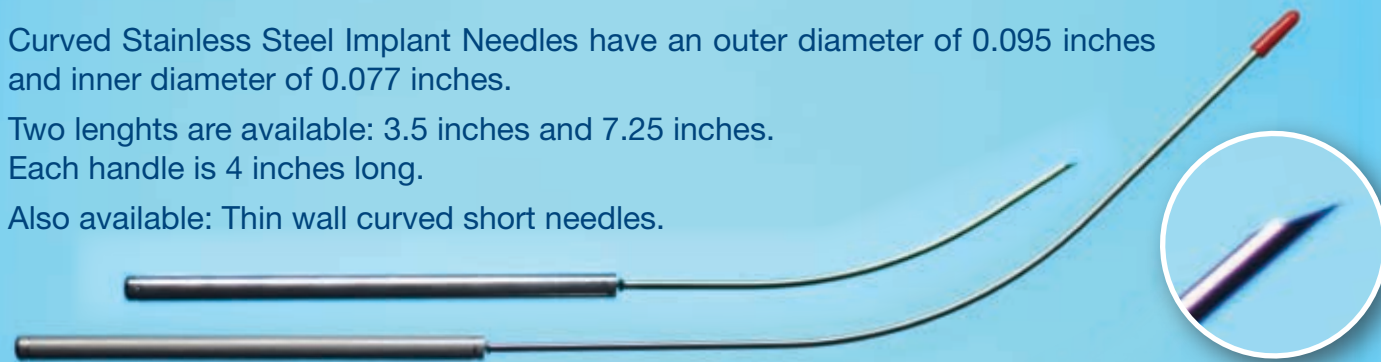

Best™ LDR Prostate Needles

Best Medical offers a wide variety of LDR Prostate Needles of different gauges, lengths and hub types. Needle kits containing customized components are available upon request.

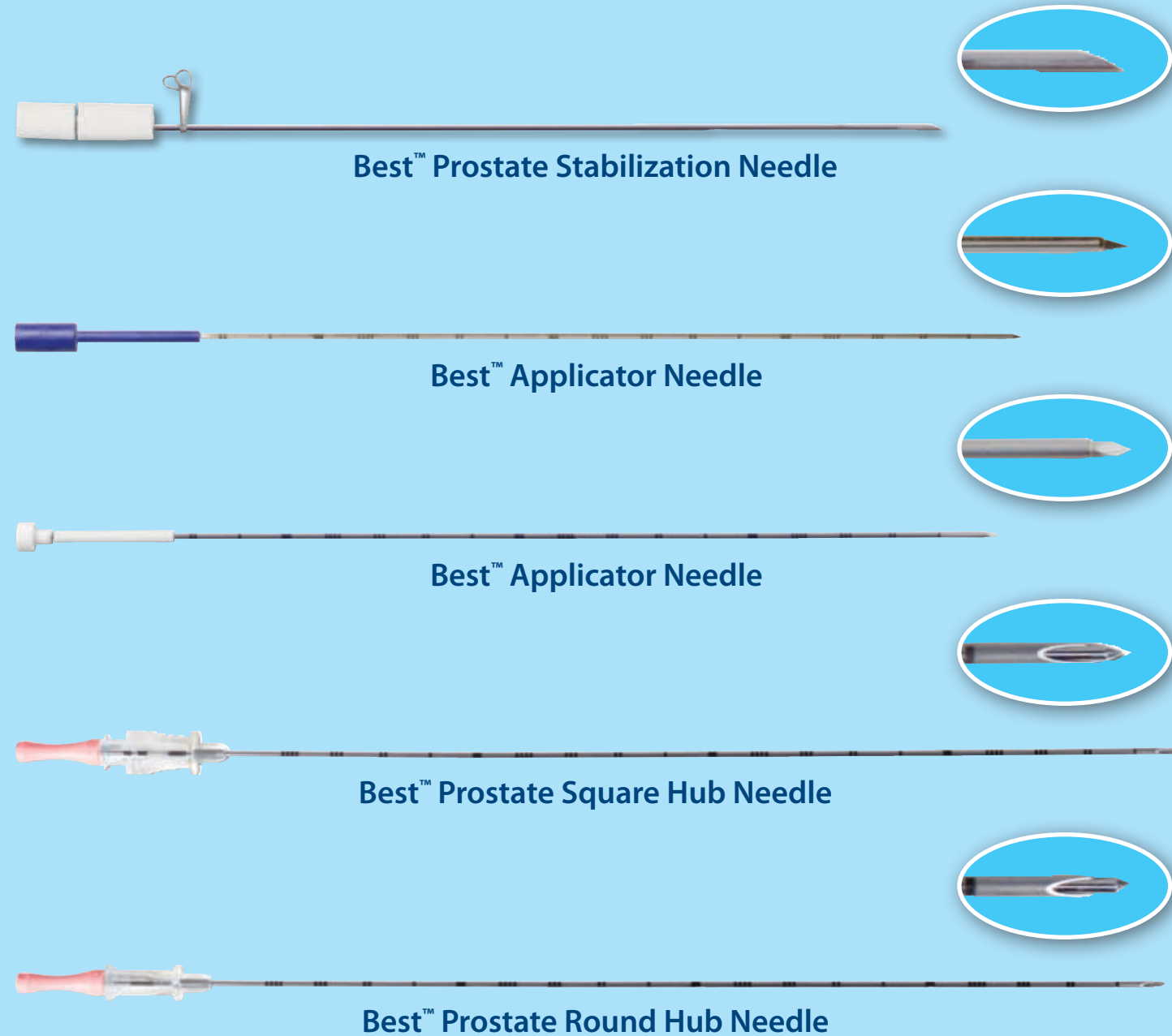

Best™ Prostate Stabilization Needle

Best™ Applicator Needle

Best™ Applicator Needle

Best™ Prostate Square Hub Needle

Best™ Prostate Round Hub Needle