


Best™ Brachytherapy Devices for GYN Applications

Best™ Localization Needles


Best provides radioactive I-125 seeds, pre-loaded in needles, sterile for localization purposes

- Best localization needles are available in 5, 7, 10 and 15 cm lengths
- Available in various seed activities with a shelf life of 90 days
- Needles are supplied in a sterile ready-to-use format and are individually packaged and shielded

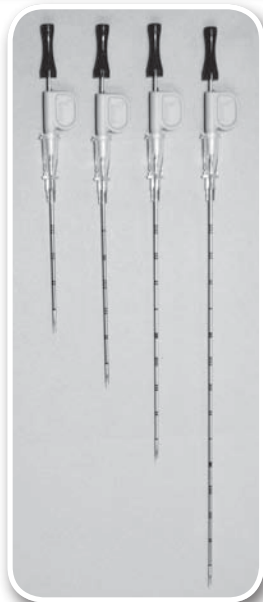

Best™ GYN Template (Disposable)

Best Medical offers a wide variety of disposable templates for HDR and LDR applications. Specify hole size when ordering (15, 16 or 17 gauge.)

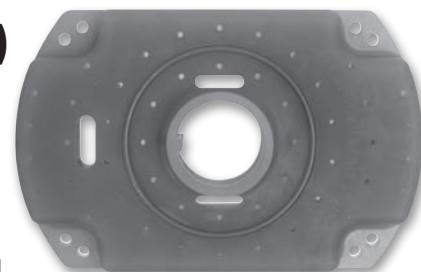

Best™ Multi-Channel Vaginal Applicator

This novel HDR applicator consists of a disposable central insert and an outer sleeve along with a positioning collar. Source channels in the central insert travel fully to the distal end to allow for optimum dose delivery.

Additionally, the outer sleeve provides a 5 mm wall, to reduce dose to the vaginal mucosa, with a 2 mm wall at the distal end of the applicator to allow for greater dose delivery. The applicator is 18 cm in length with diameters ranging from 2.0 cm to 3.5 cm.

Best™ Flexi Needles

Flexi Needles are sharp or blunt, with closed-ended needle guides, and have a nominal 0.3 cm long tip. They are made of biologically compatible material. (Custom lengths are available upon request.)



Best™ Central Rod

These cylinders are to be used with Disposable GYN and Rectal Templates. Cylinders, along with accessories, are now available with the Interstitial Disposable Vaginal Cylinder and a new "Modified" Interstitial Disposable Vaginal Cylinder. The cylinder design, when inserted into the template, positions the Flexi Needles on the surface of the cylinder.



Best Medical International, Inc. © 2020

Best Medical International, Inc. 7643 Fullerton Road, Springfield, VA 22153 USA
tel: 703 451 2378 800 336 4970 www.besttotalsolutions.com www.teambest.com

AFRICA | ASIA | EUROPE | LATIN AMERICA | MIDDLE EAST | NORTH AMERICA